KO

Lithium Disilicate Glass-ceramic Blocks & Disks

Rosetta[®] SM

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Rosetta[®] SM User Manual

Table of Contents


1. Introduction
2. Preparation Guide
3. Block Selection
4. Layer thickness
5. CAD/CAM Process
6. Before Crystallization
7. Recommended Crystallization Schedule
8. Characterizing
9. Staining & Glazing
10. Cementation
11. Workflow
12. Product Line-up
13. Indications / ContraIndications

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

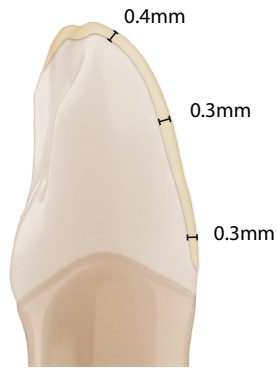
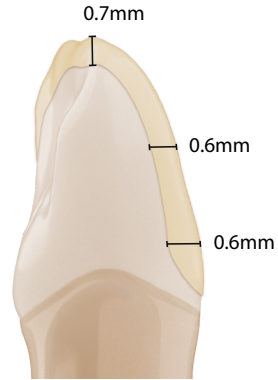
Veneer

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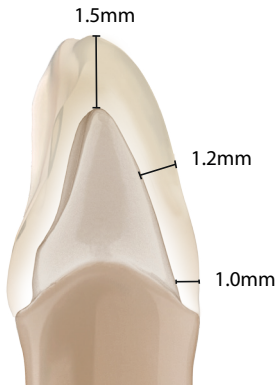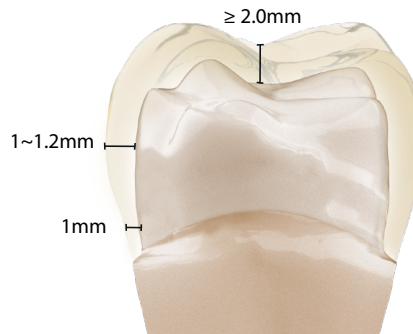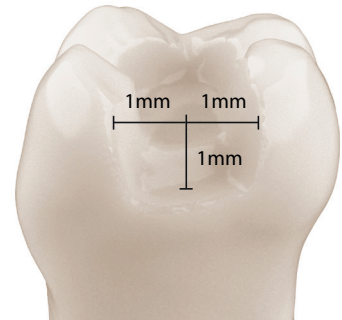
Inlay/Onlay

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Processing Technique & Indications

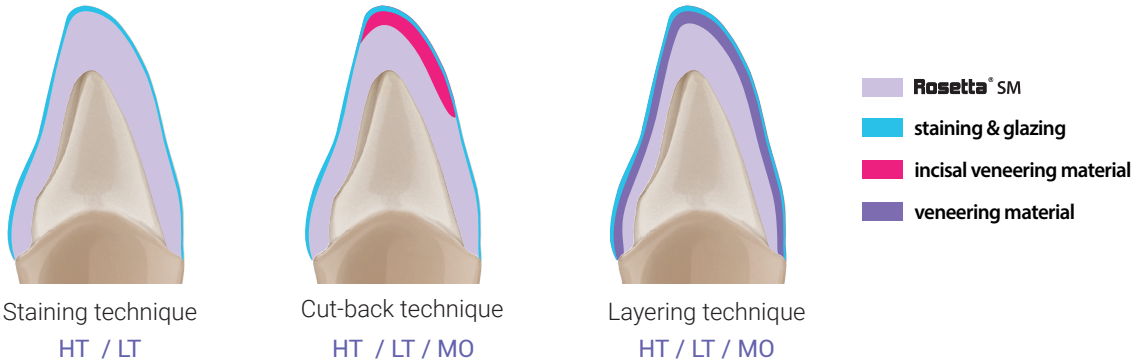
Indications								
Table Tops	Thin Veneers	Veneers	Inlays	Onlays	Partial Crowns	Anterior Crowns	Posterior Crowns	3-Unit Bridges
High Translucency								
Low Translucency								
Medium Opacity								

TIP!

For use in anterior and premolar bridge restorations can be expected.

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
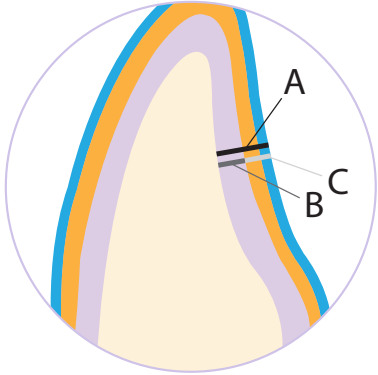
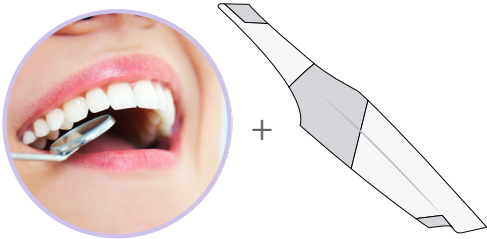
2. CAD / CAM

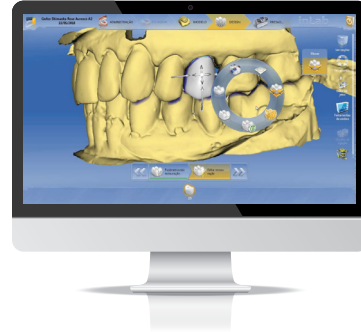
3. Wet grinding

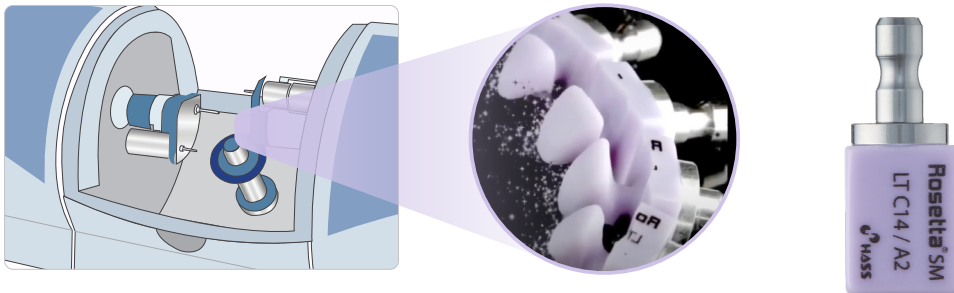
Rosetta SM is available in various CAD / CAM systems.

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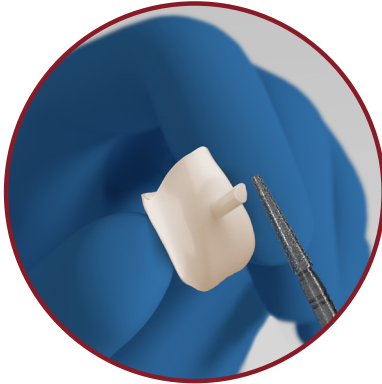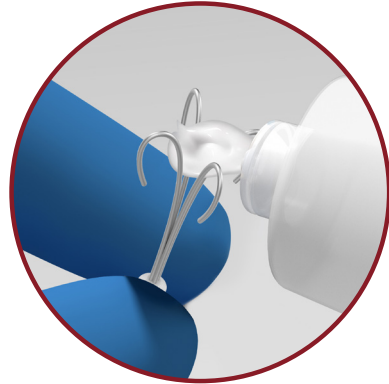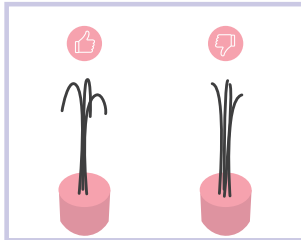
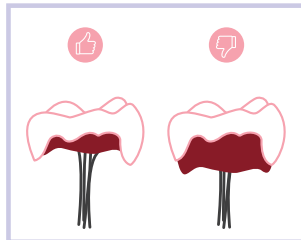
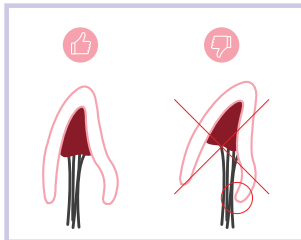
7. Recommended Crystallization Schedule

PROGRAMAT¹⁾, FOCUS²⁾

Stand-by temperature B [°C]	Closing time S [min]	Heating rate t ₁ [°C]	Firing temperature T ₁ [°C]	Holding time H ₁ [min]	VAC1 /VAC2	Long-term cooling L [°C]	Cooling time tl [°C]
400	3:00	60	840	10:00	550/840	700	0

1) PROGRAMAT is a registered trademark of IVOCCLAR VIVADENT.

2) FOCUS is a registered trademark of Shenpaz.

Please follow our recommended crystallization schedule for best results.

The name of the menu on the schedule may vary depending on the furnace's model.

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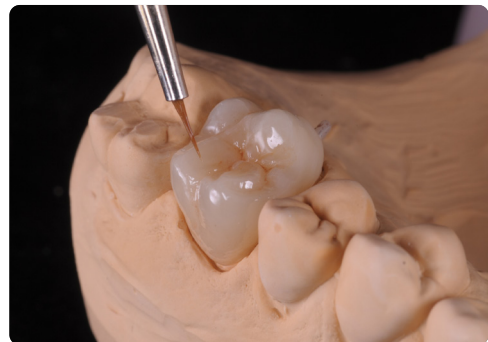
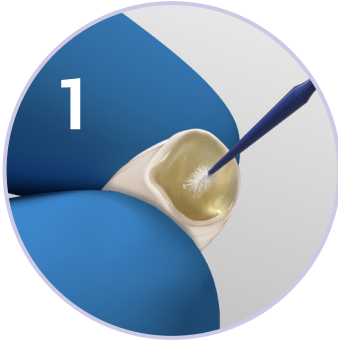
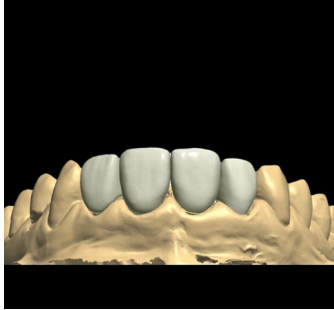
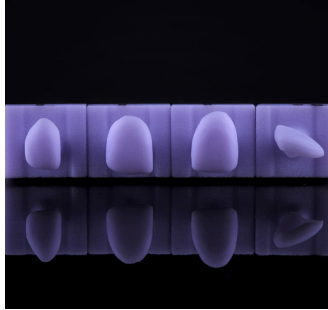
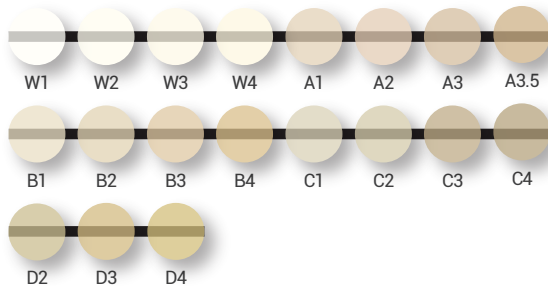
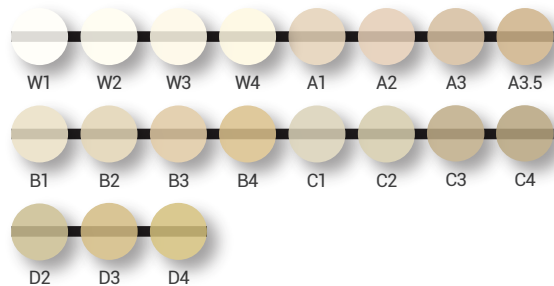
Anterior
Single Crowns

Posterior
Single Crowns

3-Unit Bridge
*up to the second Premolar

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Printed in KOREA © HASS Corporation. All rights reserved.

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